


Lipid Profiling Using the NanoMate®

Advion
BioSciences
www.advion.com
info@advion.com

This product note is based on research carried out in the laboratory of Dr. Andrej Shevchenko

Fast, reproducible, and automated
lipid analysis using chip-based
nanoelectrospray infusion

- Quantify and identify individual lipid species
- Consistent profile response from reproducible chip-based spray with up to 100 or 400 analyses per chip
- No nozzle cross contamination
- Fast and automated analysis compared to conventional nanoelectrospray

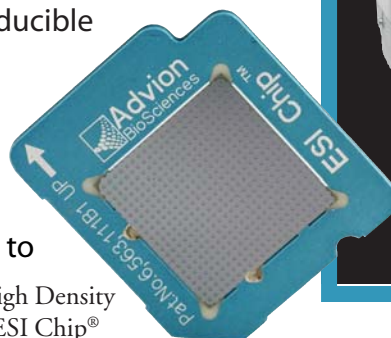


High Density
ESI Chip®



NanoMate® HD

Overlaid Precursor Ion Spectra Obtained Using NanoMate Infusion in Combination with Multiple Precursor Ion Scanning on a QSTAR Mass Spectrometer

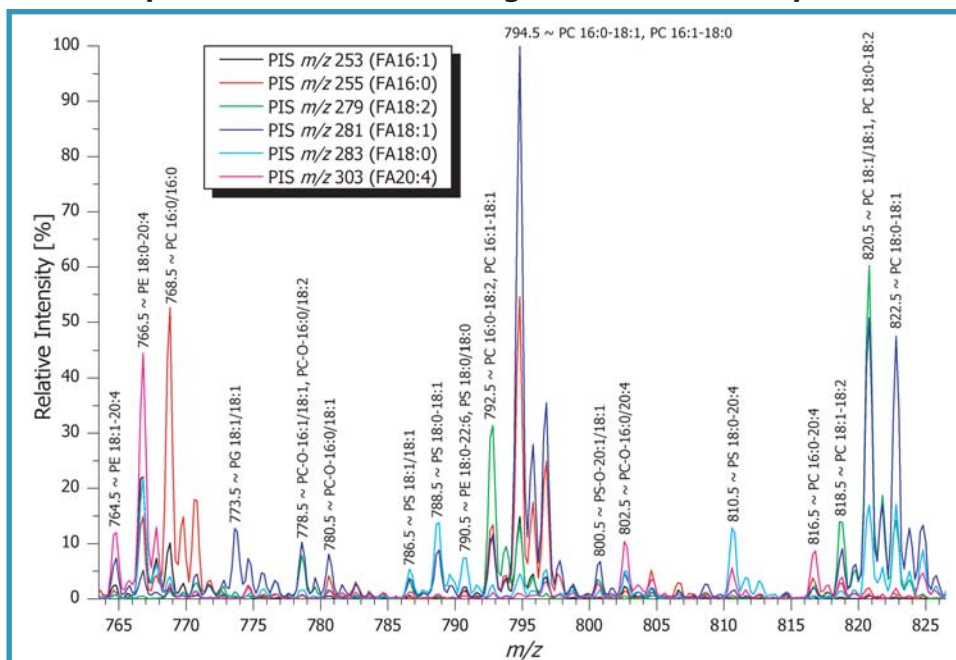

Lipid profile of mouse intestine

Glycerophospholipid species were detected by multiple precursor ion scanning for fatty acid fragment ions

Peak identification was performed using prototype Lipid Profiler software (MDS Sciex)

PC - phosphatidylcholine
PE - phosphatidylethanolamine
PS - phosphatidylserine
PG - phosphatidylglycerol
FA - fatty acid
PIS - precursor ion scan

The analysis was performed by
Christer S. Ejning

The laboratory of Dr. Andrej Schevchenko is situated at the Max-Planck Institute of Molecular Cell Biology and Genetics in Dresden, Germany. His research interests include developing novel mass spectrometry-based methodology for rapid and comprehensive quantification of individual lipid species in biological extracts.

Unleashing the Power of Mass Spectrometry™

15 Catherwood Road - Ithaca, New York - 14850
607.266.0665

26-28 Queens Road - Hethersett Norfolk - NR9 3DB
+44(0)1603 813 970